



Avoiding Entrapment

Avoiding Entrapment

Side Wedges SW



Side wedges with draw sheet. This new method of preventing a patient from falling out of bed requires NO RAILS. A draw sheet attaches to the mattress with quick release buckles. Along each side of the bed is a broad 10cm hook & loop strip. On to the strip fixes a 20cm high wedge 90cm long. The draw sheet and the wedges are covered in Vapour Permeable wipe clean fabric – Cotton drawsheet also available. The patient is able to get out of bed without the entrapment problem of rails but is also protected from the 'Roll Out'.

Also available - full length drawsheet with extra pair of wedges.

Connected Cotside Bumpers CCB



To prevent a patient from pushing limbs between the cot sides and mattress. The mattress holds the unit in place.

High Protection Bumper. To avoid entrapment and follow the requirement of the MHRA (Medicines and Healthcare products Registry Agency). Our latest cotside bumpers incorporate a highly permeable base sheet which connects the 2 cotside bumpers and forms a safety net to prevent limbs falling through under the rails.

CCB 140 x 76 cm
CCBPB 200 x 76 cm
CCBPH 200 x 90cm

Fallout Mat FOM



Heavy duty dense foam crash mat can be used with standard bed to low level beds. Folds easily for storing – cover is non slip wipe clean material.

Also available with straps at each end.
Size 198 x 60 x 5 cm

Cotside Bumpers CSB (Guardian Angels)




Designs suitable for use in children's units

To fit a vast range of bed rails from 140cm to 205cm long and various heights. The Basic wipe clean bumper has a foam inner and full length hook and loop fastenings along the bottom edge. Available in extra high to fit high rails when an extra overlay is being used to fully enclose extending trombone rails.



Quilted Hollow Fibre Filled Cotside Bumpers CSS, CSL

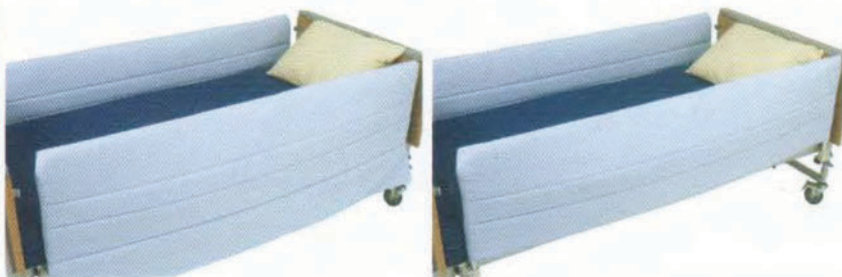

Covered in washable F.R. polyester, with extra long tucks to secure under mattress to prevent limbs falling between the rail and mattress.

CSS 140 x 43 cm or 54 x 17 inches
CSL 183 x 43 cm or 72 x 17 inches

Quilted - Wipe Clean Hollow Fibre Filled Bumpers CSQS, CSQL



Made with closed ends or open ended in various lengths these are available in fully washable and wipe clean material - see price list for details.

CSQS 140 x 43 cm
CSQL 183 x 43 cm

Deep Pile Fleece Bumpers



When the patient is subjected to fits – this model will allow air to pass through, see price list for details.

CSHP1 137 x 43 cm or 54 x 17 inches
CSHP2 173 x 43 cm or 54 x 17 inches

See Through Side Rail Protectors (Connected)




- Enables the patient to look around
- Removes any feeling of claustrophobia
- More importantly! Clear sides allow continuous patient monitoring from a distance.

Specification: A rounded foam pad covered in wipe clean material 'clips' onto the top rail and a strong clear PVC panel closes the spaces between the side rails. The protectors are held firmly in place by 3 or 4 vertical hook and loop straps, depending on the length: 135cm - 183cm - 200cm.

Extra Secure Model has a net which passes under the mattress and connects the two side panels together, thus closing the gap between the rails and the bed frame.

Net to Close Gaps in Rails NET

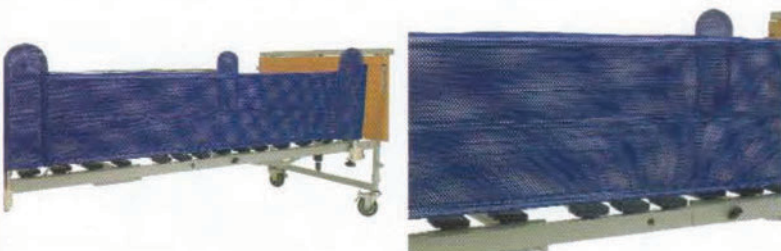

A simple strong net covers the rails – tea cosy style – tough and closed along the bottom edge.